

Silicone Super-Absorbent Fiber Dressings



Vitale® Silicone Super-Absorbent Dressings are highly absorbent dressings for wounds with light to heavy drainage. The absorbent pad is surrounded by a waterproof but vapor permeable border. The adhesive is soft, flexible silicone for maximum patient comfort & nursing convenience.



This dressing is composed of 4 important layers:

- 1) The Silicone Contact Layer covers the entire patient contact surface. It allows transmission of fluids to the absorbent pad, while maintaining the ideal moist environment for healing. Silicone reduces pain and trauma during dressing application, inspection, and at change time.
- 2) The perforated Silicone film behind the Silicone Adhesive Contact Layer provides firm support for the silicone adhesive and has many holes in it to allow moisture to penetrate.
- 3) The highly absorbent super-absorbent fiber pad rapidly absorbs exudates and holds them until dressing change time.
- 4) The outer film, breathable yet waterproof, is a barrier to bacteria and other microbes.

These dressings come in four sizes plus 2 sizes of sacral dressings to accommodate all sizes of wounds.

Vitale Silicone Super-Absorbent Dressings may be used as a primary dressing in direct contact with the wound bed, or as a secondary (cover) dressing over an appropriate wound filler or other primary dressing.


Indications: Vitale Silicone Super-Absorbent Dressings are for single use only and can be used for the treatment of the following wounds:

- | | | |
|---|--------------------|----------------------------------|
|  Superficial to Deep | • Pressure Ulcers | • Second Degree Burns |
|  Light to Heavy Exuding Wounds | • Diabetic Ulcers | • Plus Covering Other Treatments |
| | • Arterial Ulcers | |
| | • Skin donor sites | |

The dressings may remain in place for 1-7 days based on wound characteristics and drainage amounts. Not indicated for dry wound conditions or for 3rd degree burns. Vitale Silicone Super-Absorbent Dressings are sterilized by exposure to ethylene oxide gas and must not be resterilized.



Two sacral sizes available.

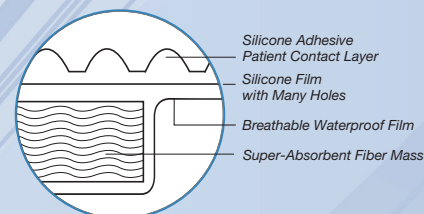

Vitale Silicone Super-Absorbent Dressings use a unique release material layout, with three sections instead of the usual two. This allows the user to take off the center release paper first, and position the dressing properly using the silicone adhesive to hold the dressing in place. Then the side release material can then be easily removed to complete application.

– Ordering Information –

Cat. #:	Description:	Size:	Qty/box:	HCPC Code:
20234	Vitale Silicone Super-Absorbent Dressings	3½" x 4" (9x10 cm)	10	A6196
20245	Vitale Silicone Super-Absorbent Dressings	4¾" x 5½" (12x14 cm)	10	A6196
20267	Vitale Silicone Super-Absorbent Dressings	6" x 7" (15x18 cm)	10	A6197
20280	Vitale Silicone Super-Absorbent Dressings	8" x 10" (20x25 cm)	10	A6197
20290	Vitale Sacral Silicone Super-Absorbent Dressings, Large	7" x 6¾" (175x170 mm)	5	A6197
20292	Vitale Sacral Silicone Super-Absorbent Dressings, Small	6" x 5.8" (15x14.5 cm)	10	A6196

Features & Benefits:

Highly absorbent
Multi-functional
Moist wound environment
Silicone adhesive
Thin profile

Absorbs many times its own weight, allowing the dressing to remain in the wound.
May be used as primary or secondary dressing as appropriate.
The 4-layer pad keeps the wound moist for healing.
Won't irritate peri-wound tissue. Border can be lifted to inspect wound site between dressing changes.
More comfortable; less stiff.

