

VITÁLE®

Silicone Foam Dressings




Vitále® Silicone Foam Dressings are highly absorbent dressings for wounds with moderate to heavy drainage. The absorbent pad is surrounded by a vapor permeable border. The adhesive is soft, flexible silicone for maximum patient comfort & nursing convenience.

This dressing is composed of 5 important layers:

- 1) The Silicone Contact Layer covers the entire patient contact surface. It allows transmission of fluids to the foam, while helping maintain the ideal moist environment for healing. Silicone reduces pain and trauma during dressing application, inspection, and at change time.
- 2) The perforated Silicone film behind the Silicone Adhesive Contact Layer provides firm support for the adhesive and allows moisture to penetrate.
- 3) The highly absorbent foam pad rapidly absorbs exudates and holds them until dressing change time.
- 4) Behind the foam is a Super-Absorbent Fiber (SAF) layer that pulls exudates into the dressing and keeps them from soaking back into the wound site.

5) The outer film, breathable yet waterproof, is a barrier to bacteria and viruses.

Vitále Foam Dressings are ideal for draining wounds, reducing the risk of maceration. The 5-layer pad construction keeps new granulation tissue moist, prevents it from drying out, and encourages epithelialization.

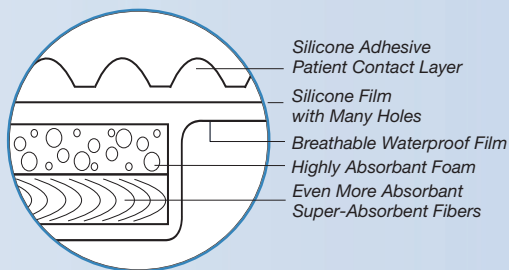

Vitále Sacral Silicone Foam Dressings are now available...a convenient shape for ease of application & comfort in the sacral area.



Indications: Vitále Silicone Foam Dressings are for single use only and can be used for the treatment of the following wounds:

- Pressure Ulcers
- Arterial Ulcers
- Venous Ulcers
- Skin donor sites
- Superficial & Deep Wounds
- Moderate to Heavy Exudating Wounds
- Diabetic Ulcers
- Second Degree Burns
- Plus a Special Sacral Size

Not indicated for dry wound conditions or for 3rd degree burns. Vitále Silicone Foam Dressings are sterilized by exposure to ethylene oxide gas and must not be resterilized.

– Ordering Information –

Cat. #:	Description:	Size:	Qty/box:	HCPC Code:
20133	Vitále Silicone Foam Dressings	3"x3" (7.5x7.5 cm)	10	A6212
20144	Vitále Silicone Foam Dressings	4"x4" (10x10 cm)	10	A6212
20166	Vitále Silicone Foam Dressings	6"x6" (15x15 cm)	10	A6212
20190	Vitále Sacral Silicone Foam Dressings	7"x6.75" (17.5x17cm)	5	A6212

Features & Benefits:

- Highly absorbent *Absorbs many times its own weight, allowing the dressing to remain in the wound for up to 7 days.*
- Moist wound environment *The 5-layer pad keeps the wound moist for healing.*
- Silicone adhesive *Sticks firmly but won't irritate peri-wound tissue. Can be lifted to inspect between dressing changes.*



CellEra LLC

Cincinnati, OH

www.cellerallc.com

800-892-2391